


MultiHarp 150

High-Throughput Multichannel Event Timer & TCSPC Unit

- 4 or 8 independent input channels
- one common sync channel (up to 1.2 GHz)
- High sustained data throughput (80 Mcps in time tagging mode, 180 Mcps in histogramming mode)
- 80ps base resolution
- Ultra short dead time (650 ps) per channel, no dead time across channels
- USB 3.0 interface
- Marker inputs for external synchronization signals or other events
- White Rabbit ready

Applications

- Time-resolved fluorescence and luminescence spectroscopy
- Coincidence correlation, antibunching
- Quantum optics
- FLIM, FRET, FCS
- Multicolor lifetime imaging
- Time-of-Flight (ToF) measurements, LIDAR
- Single Molecule Spectroscopy (SMS)

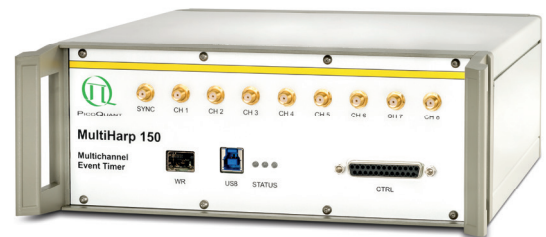

The MultiHarp 150 is an easy-to-use, plug-and-play multi-channel event timer and Time-Correlated Single Photon Counting (TCSPC) device. Its outstanding features are an extremely fast signal processing along with an extraordinarily high sustained data throughput via USB 3.0 interface. The MultiHarp 150 is a compact, robust and reliable device whose high quality is reflected by our unique 5-year limited warranty.

Multiple input channels, outstanding flexibility

Two MultiHarp 150 models are available with either 4 or 8 identical detection channels, which are synchronized but independent. Each model features also one common synchronization input. All channels including the sync input can be used as detector inputs, e.g. for coincidence correlation or coincidence counting. The MultiHarp 150 is also ideally suited for performing TCSPC with multiple detectors using forward start-stop operation. Here the common sync channel allows for synchronization with the excitation source.

Ultrashort dead time for high data throughput

The smartly designed timing electronics allows to fully exploit the count rate limits of time-correlated single photon counting, without having to compromise on the time resolution for many modern single photon detectors. The ultra short dead time of 650 ps allows detecting multiple photons per excitation cycle even at the highest repetition rates achievable by modern picosecond pulsed lasers (requires a detector from the PMA Hybrid Series).

Adjustable timing offsets for each input channel

Each input channel has internal adjustable timing offsets with ± 100 ns range at 80 ps steps. This greatly simplifies the cabling requirements for experimental set-ups, as selecting cable lengths to compensate for signal delays is no longer needed.

Time-Tagged Time-Resolved modes


The Time-Tagged Time-Resolved (TTTR) modes supported by the MultiHarp 150 record all relevant time and channel routing information of each detected individual photon event. By storing this full data set, it becomes possible to carry out the most comprehensive and sophisticated analysis of photon dynamics after a measurement. A real time data correlator is included to monitor FCS experiments at count rates of up to 1 000 000 counts/sec. Furthermore, the MultiHarp 150 can be synchronized with other hardware such as scanners when operated in TTTR mode.

White Rabbit ready

White Rabbit is a fully deterministic, Ethernet-based timing network which provides sub-nanosecond accuracy and precise synchronization of devices over large distances. Thanks to its White Rabbit interface, the MultiHarp 150 is prepared to be used in networks that are based on this emerging technology.

Easy to use software included, custom programming supported

The MultiHarp 150 comes with a Windows software package that providing all important functions such as setting measurement parameters, displaying results, loading / saving of measurement parameters and measurement curves. Important measurement data, including count rate, count maximum, position and peak width are continuously displayed. A comprehensive online help system eases the user into fully employing the capabilities of the MultiHarp 150. A library for custom programming, e.g., with LabVIEW is also included. PicoQuant is committed to the support and development of this software and upgrades with extended functionality will be made available.

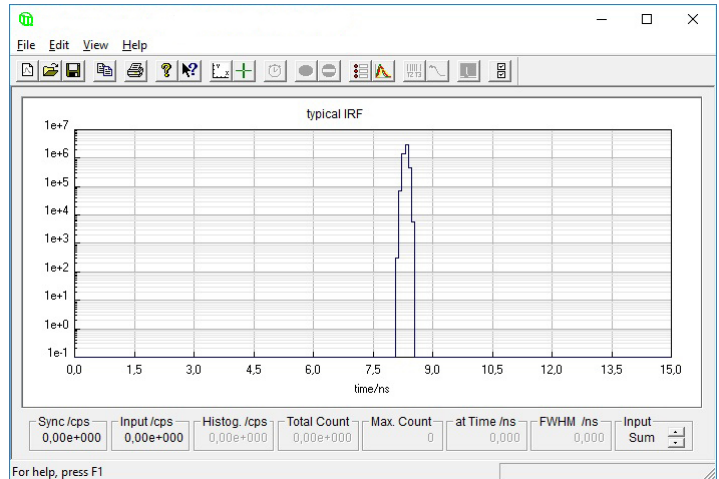

The MultiHarp 150 can be used in different operation modes:

Integration mode, oscilloscope mode, Time-Tagged Time-Resolved mode (listing each event arrival time), on-line („real-time“) correlator for FCS

Specifications

Input Channels and Sync	Constant level trigger on all inputs, software adjustable
Number of detector channels	4 (MultiHarp 150 4N) or 8 (MultiHarp 150 8N) plus sync
Input voltage operating range (pulse peak into 50 Ohms)	- 1200 mV to 1200 mV
Input voltage max. range (damage level)	± 2500 mV
Trigger edge	falling or rising edge, software adjustable
Trigger pulse width	> 0.4 ns (rise time max. 20 ns)
Time to Digital Converters	
Time bin width (adjustable)	80 ps, 160 ps, 320 ps, [...], 335.5 µs
Electrical time resolution	< 85 ps rms
Dead time	< 650 ps
Adjustable programmable time offset for each input channel	± 100 ns, resolution 80 ps
Differential non-linearity	< 15 % peak, < 1.5 % rms (over full measurement range)
Maximum sync rate (periodic pulse train)	1.2 GHz

Histogrammer	
Count depth	32 bit (4 294 967 296 counts)
Full scale time range	5.24 μ s to 21.99 s
Maximum number of time bins	65536
Peak count rate per input channel	1.5×10^9 counts/sec for 1 200 events
Total sustained count rate, sum over all input channels	180×10^6 counts/sec
TTR Engine	
T2 mode resolution	80 ps
T3 mode resolution	80 ps, 160 ps, 320 ps, [...], 335.5 μ s
FiFo buffer depth (records)	134 217 728 events
Peak count rate per input channel	1.5×10^9 counts/sec for 1 200 events
Total sustained count rate, sum over all input channels	80×10^6 counts/sec
Operation	
PC interface	USB 3.0
PC requirements	Dual core CPU or better, min. 2 GHz CPU clock, min. 4 GB memory
Operating system	Windows 8/10
Power consumption	50 W



PicoQuant GmbH
Rudower Chaussee 29 (IGZ)
12489 Berlin
Germany

Phone +49-(0)30-1208820-0
Telefax +49-(0)30-1208820-90
Email info@picoquant.com
Web www.picoquant.com